

Evaporation and Intermolecular Attractions

In this experiment, Temperature Probes are placed in various liquids. Evaporation occurs when the probe is removed from the liquid's container. This evaporation is an endothermic process that results in a temperature decrease. The magnitude of a temperature decrease is, like viscosity and boiling temperature, related to the strength of intermolecular forces of attraction. In this experiment, you will study temperature changes caused by the evaporation of several liquids and relate the temperature changes to the strength of intermolecular forces of attraction. You will use the results to predict, and then measure, the temperature change for several other liquids.

You will encounter two types of organic compounds in this experiment—alkanes and alcohols. The two alkanes are pentane, C_5H_{12} , and hexane, C_6H_{14} . In addition to carbon and hydrogen atoms, alcohols also contain the $-OH$ functional group. Methanol, CH_3OH , and ethanol, C_2H_5OH , are two of the alcohols that we will use in this experiment. You will examine the molecular structure of alkanes and alcohols for the presence and relative strength of two intermolecular forces—hydrogen bonding and dispersion forces.

OBJECTIVES

In this experiment, you will

- Study temperature changes caused by the evaporation of several liquids.
- Relate the temperature changes to the strength of intermolecular forces of attraction.

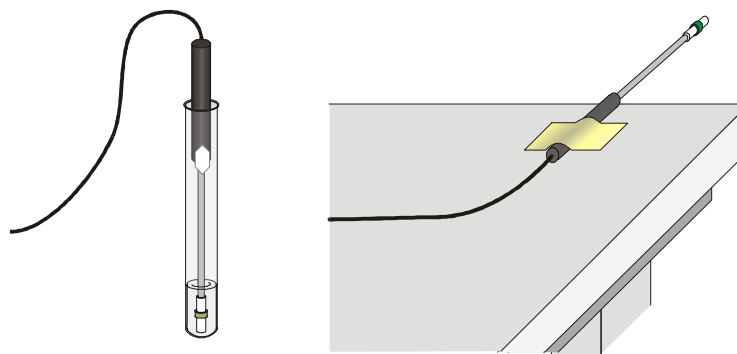


Figure 1

MATERIALS

LabPro interface	methanol (methyl alcohol)
Palm handheld	ethanol (ethyl alcohol)
Data Pro program	1-propanol
2 Temperature Probes	1-butanol
6 pieces of filter paper (2.5 cm x 2.5 cm)	n-pentane
2 small rubber bands	n-hexane
masking tape	

PRE-LAB EXERCISE

Prior to doing the experiment, complete the Pre-Lab table. The name and formula are given for each compound. Draw a structural formula for a molecule of each compound. Then determine the molecular weight of each of the molecules. Dispersion forces exist between any two molecules, and generally increase as the molecular weight of the molecule increases. Next, examine each molecule for the presence of hydrogen bonding. Before hydrogen bonding can occur, a hydrogen atom must be bonded directly to an N, O, or F atom within the molecule. Tell whether or not each molecule has hydrogen-bonding capability.

PROCEDURE

1. Obtain and wear goggles! **CAUTION:** *The compounds used in this experiment are flammable and poisonous. Avoid inhaling their vapors. Avoid contacting them with your skin or clothing. Be sure there are no open flames in the lab during this experiment. Notify your teacher immediately if an accident occurs.*
2. Plug Temperature Probe 1 into Channel 1 and Temperature Probe 2 into Channel 2 of the LabPro interface. Connect the handheld to the LabPro using the interface cable. Firmly press in the cable ends.
3. Turn on the handheld. To start Data Pro, tap the Data Pro icon on the Applications screen. Choose New from the Data Pro menu or tap **[New]** to reset the program.
4. Set up the data-collection mode for two Temperature Probes.
 - a. Tap **[Setup]** on the Main screen.
 - b. If the handheld displays TEMP(C) in CH 1 and CH 2, proceed directly to Step 5. If it does not, continue with this step to set up your sensors manually.
 - c. Tap **[CH1:]** to select Channel 1.
 - d. Press the Scroll buttons on the handheld to scroll through the list of sensors.
 - e. Choose the Temperature Probe you are using (in °C) from the list of sensors.
 - f. Tap **[CH2:]** to select Channel 2.
 - g. Press the Scroll buttons on the handheld to scroll through the list of sensors.
 - h. Choose the Temperature Probe you are using (in °C) from the list of sensors.
5. Set up the data-collection mode for two Temperature Probes.
 - a. On the Setup screen, tap **[Settings:]**.
 - b. Enter “3” as the time between samples in seconds, using the onscreen keyboard (tap “123”) or using the Graffiti writing area.
 - c. Enter “80” as the number of samples. (Data collection will last 4 minutes.)
 - d. Tap **[OK]** twice to return to the Main screen.
6. Wrap Probe 1 and Probe 2 with square pieces of filter paper secured by small rubber bands as shown in Figure 1. Roll the filter paper around the probe tip in the shape of a cylinder. Hint: First slip the rubber band on the probe, wrap the paper around the probe, and then finally slip the rubber band over the paper. The paper should be even with the probe end.
7. Stand Probe 1 in the ethanol container and Probe 2 in the methanol container. Make sure the containers do not tip over.

8. Prepare 2 pieces of masking tape, each about 10 cm long, to be used to tape the probes in position during Step 9.
9. After the probes have been in the liquids for at least 30 seconds, tap to begin collecting temperature data. A live graph of temperature vs. time for both Probe 1 and Probe 2 is being plotted on the handheld screen. The live readings are displayed in the upper-right corner of the graph, Probe 1 first, Probe 2 below. Monitor the temperature for 15 seconds to establish the initial temperature of each liquid. Then simultaneously remove the probes from the liquids and tape them so the probe tips extend 5 cm over the edge of the table top as shown in Figure 1.
10. Data collection will stop after 4 minutes (or tap if you want to end data collection before 4 minutes has elapsed). On the displayed graph of temperature vs. time, each point for Probe 1 is plotted with open squares, and each point for Probe 2 with filled squares.
11. To examine the data pairs on the displayed graph, tap or any data point. As you move the examine line, the temperature values of both Probe 1 (top value) and Probe 2 (lower value) are displayed to the right of the graph. Based on your data, determine the maximum temperature, t_1 , and minimum temperature, t_2 . Record t_1 and t_2 for each probe.
12. For each liquid, subtract the minimum temperature from the maximum temperature to determine Δt , the temperature change during evaporation.
13. Based on the Δt values you obtained for these two substances, plus information in the Pre-Lab exercise, *predict* the size of the Δt value for 1-butanol. Compare its hydrogen-bonding capability and molecular weight to those of ethanol and 1-propanol. Record your predicted Δt , then explain how you arrived at this answer in the space provided. Do the same for n-pentane. It is not important that you predict the exact Δt value; simply estimate a logical value that is higher, lower, or between the previous Δt values.
14. Test your prediction in Step 13 by repeating Steps 7-12 using 1-butanol with Probe 1 and n-pentane with Probe 2.
15. Based on the Δt values you have obtained for all four substances, plus information in the Pre-Lab exercise, predict the Δt values for methanol and n-hexane. Compare the hydrogen-bonding capability and molecular weight of methanol and n-hexane to those of the previous four liquids. Record your predicted Δt , then explain how you arrived at this answer in the space provided.
16. Test your prediction in Step 15 by repeating Steps 7-12, using methanol with Probe 1 and n-hexane with Probe 2.

PROCESSING THE DATA

1. Two of the liquids, n-pentane and 1-butanol, had nearly the same molecular weights, but significantly different Δt values. Explain the difference in Δt values of these substances, based on their intermolecular forces.
2. Which of the alcohols studied has the strongest intermolecular forces of attraction? The weakest intermolecular forces? Explain using the results of this experiment.
3. Which of the alkanes studied has the stronger intermolecular forces of attraction? The weaker intermolecular forces? Explain using the results of this experiment.

Experiment 9

4. Plot a graph of Δt values of the four alcohols versus their respective molecular weights. Plot molecular weight on the horizontal axis and Δt on the vertical axis.

PRE-LAB

Substance	Formula	Structural Formulas	Molecular Weight	Hydrogen Bond (Yes or No)
ethanol	C ₂ H ₅ OH			
1-propanol	C ₃ H ₇ OH			
1-butanol	C ₄ H ₉ OH			
n-pentane	C ₅ H ₁₂			
methanol	CH ₃ OH			
n-hexane	C ₆ H ₁₄			

DATA TABLE

Substance	t ₁ (°C)	t ₂ (°C)	Δt (t ₁ -t ₂) (°C)	Predicted Δt (°C)	Explanation
ethanol					
1-propanol					
1-butanol					
n-pentane					
methanol					
n-hexane					

TEACHER INFORMATION

Evaporation and Intermolecular Attractions

1. We recommend wrapping the probes with paper as described in the procedure. Wrapped probes provide more uniform liquid amounts, and generally greater Δt values, than bare probes. Chromatography paper, filter paper, and various other paper types work well.
2. Snug-fitting rubber bands can be made by cutting short sections from a small rubber hose. Surgical tubing works well. Orthodontist's rubber bands are also a good size.
3. Other liquids can be substituted. Although it has a somewhat larger Δt , 2-propanol can be substituted for 1-propanol. Some petroleum ethers have a high percentage of hexane and can be used in its place. Other alkanes of relatively high purity, such as n-heptane or n-octane can be used. Water, with a Δt value of about 5°C, emphasizes the effect of hydrogen bonding on a low-molecular weight liquid. However, students might have difficulty comparing its hydrogen bonding capability with that of the alcohols used.
4. Sets of the liquids can be supplied in 13 × 100-mm test tubes stationed in stable test-tube racks. This method uses very small amounts of the liquids. Alternatively, the liquids can be supplied in sets of small bottles kept for future use. Adjust the level of the liquids in the containers so it will be above the top edge of the filter paper.
5. Because several of these liquids are highly volatile, keep the room well-ventilated. Cap the test tubes or bottles at times when the experiment is not being performed. The experiment should not be performed near any open flames.
6. Other properties, besides Δt values, vary with molecular size and consequent size of intermolecular forces of attraction. Viscosity increases noticeably from methanol through 1-butanol. The boiling temperatures of methanol, ethanol, 1-propanol, and 1-butanol are 65°C, 78°C, 97°C, and 117°C, respectively.
7. If you have only one Temperature Probe per interface, this experiment can still be completed in one class period. It is also possible to do four of the liquids during one class period, and the remaining two liquids the next day. This provides students with additional time to consider their predictions. The Stainless Steel Temperature Probe and the Direct-Connect Temperature Probe work very well for this experiment. We do not recommend that you use the TI Temperature Probe that was shipped with the original CBL; our tests show that pentane liquid sometimes penetrates the seal on the tip of the TI-Temperature Probe.
8. HAZARD ALERTS:
 - n-Hexane: Flammable liquid; dangerous fire risk; may be irritating to respiratory tract. Hazard Code: B—Hazardous.
 - Methanol: Flammable; dangerous fire risk; toxic by ingestion (ingestion may cause blindness). Hazard Code: B—Hazardous.
 - Ethanol: Dangerous fire risk; flammable; addition of denaturant makes the product poisonous—it cannot be made non-poisonous; store in a dedicated flammables cabinet or safety cans. If a flammables cabinet or safety cans are not available, store in a Flinn *Saf-Stor*® Can. Hazard Code: C—Somewhat hazardous.

Experiment 9

n-Pentane: Flammable liquid; narcotic in high concentrations. Hazard Code: B—Hazardous.

1-Propanol: Flammable liquid; dangerous fire risk; harmful to eyes and respiratory tract. Hazard Code: B—Hazardous.

1-Butanol: Moderate fire risk; toxic on prolonged inhalation; eye irritant; absorbed by skin. Substance may develop explosive hydroperoxides. Use fresh materials only. All Flinn chemicals are date labeled. Hazard Code: B—Hazardous.

The hazard information reference is: Flinn Scientific, Inc., *Chemical & Biological Catalog Reference Manual*, 2000, P.O. Box 219, Batavia, IL 60510. See *Appendix F* of this book, *Chemistry with Calculators*, for more information.

ANSWERS TO QUESTIONS

Answers have been removed from the online versions of Vernier curriculum material in order to prevent inappropriate student use. Graphs and data tables have also been obscured. Full answers and sample data are available in the print versions of these labs.

DATA TABLE


Substance	t_1 (°C)	t_2 (°C)	Δt ($t_1 - t_2$) (°C)
ethanol	xxxx	xxxx	xxxx
1-propanol	xxxx	xxxx	xxxx
1-butanol	xxxx	xxxx	xxxx
n-pentane	xxxx	xxxx	xxxx
Methanol	xxxx	xxxx	xxxx
n-hexane	xxxx	xxxx	xxxx

Predicted Δt (°C)	Explanation
varies (< 4.9°C)	xxxx
varies (> 8.3°C)	xxxx
varies (> 8.3°C)	xxxx
varies (< 16.1°C)	xxxx


PRE-LAB RESULTS

Substance	Formula	Structural Formulas	Molecular Weight	Hydrogen Bond (Yes or No)
ethanol	XXXX	<pre> X X X-X-X-X-X X X </pre>	XXXX	XXXX
1-propanol	XXXX	<pre> X X X-X-X-X-X X X </pre>	XXXX	XXXX
1-butanol	XXXX	<pre> X X X-X-X-X-X X X </pre>	XXXX	XXXX
n-pentane	XXXX	<pre> X X X-X-X-X-X X X </pre>	XXXX	XXXX
methanol	XXXX	<pre> X X X-X-X-X-X X X </pre>	XXXX	XXXX
n-hexane	XXXX	<pre> X X X-X-X-X-X X X </pre>	XXXX	XXXX


Experiment 9

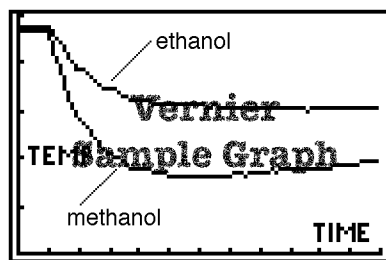


Evaporation of methanol and ethanol

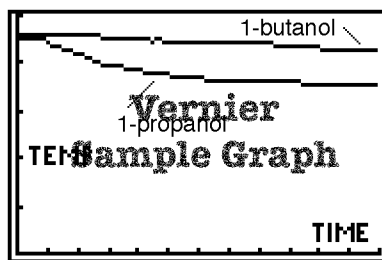

Evaporation of 1-propanol and 1-butanol

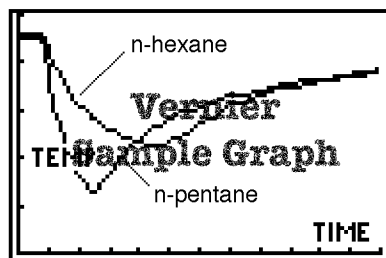

Evaporation of n-pentane and n-hexane

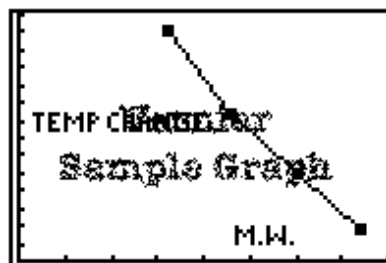

Temperature change vs. alcohol molecular wt.